
MaxiGo 2AE

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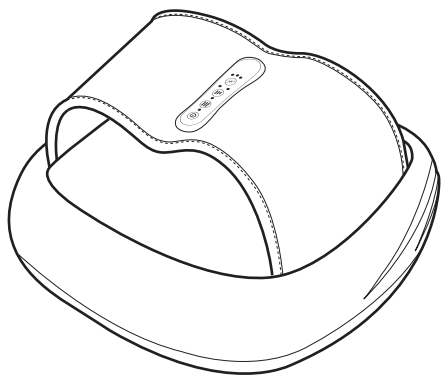Адрес производителя: D Building No. 100, JINFU ROAD, Chengnan Industrial Zone, TONG'AN District, Xiamen, China

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
foot massager

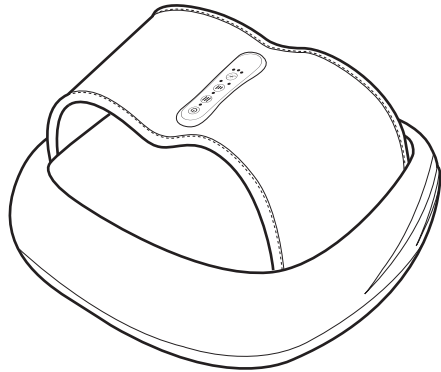Thank you for

- Choosing this product to complement your healthcare needs.
- Before using please read this user manual carefully for the correct usage and operation of this product.
- Please keep the user manual on hand for easy reference.

Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER. To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING. To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 4) Do not drop or insert any object into any opening.
- 5) Do not use outdoors.
- 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
- 8) Do not carry this product by its supply cord or use cord as a handle.
- 9) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
- 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.

11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination.

12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.

13) Use this appliance only for its intended use as described in this manual.

14) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

15) The appliance is only to be used with the power supply unit provided with the appliance.

16) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

- Children should be supervised to ensure that they do not play with the appliance.

- Cleaning and user maintenance shall not be made by children without supervision.

17) As the surface of this product may get hot, people who are sensitive to heat should pay attention to using it to avoid minor burns.

- If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health. 1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.

2. People who are supposed to lose their sense of temperature.

3. People with malignant tumors.

4. People with heart disease.

5. A person who is restless early in pregnancy or shortly after giving birth.

6. One who has an abnormal or curved spine.

7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.

8. Someone with a scar on their skin.

9. People who need quiet and people who are in poor health.

10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.

11. People with 38 °C temperature (fever).

12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.

13. People with tendonitis.

14. Others treated at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)



• Do not wash!

• Do not dry clean!

Name and Function of The Machine Parts

Configuration Figure of Functions

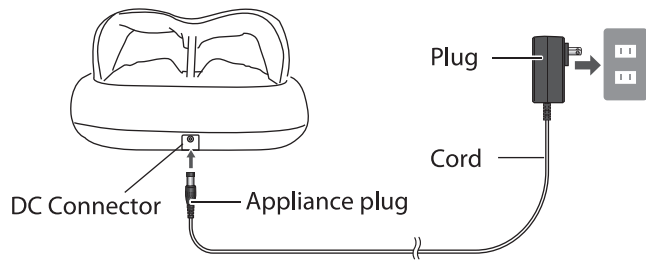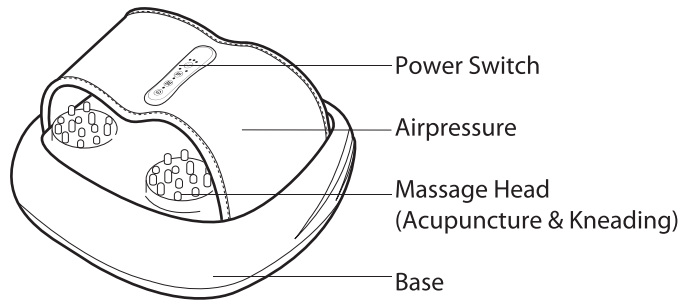How to Operate

1. Power ON/OFF:

Press this button, machine power-on.

All keys are effective. Automatically setting time is 15 mins.

The Acupuncture and Air pressure function turn on when power on, and the air-pressure of the feet is opened by default in low intensity, acupuncture & kneading speed defaults to slow speed.

Press this button again, the whole machine is shut down, and all massages are stopped and the system is in standby mode.

When time's up, machine stop massaging, enter standby state.

2. Speed:

It is defaulted to be in Slow speed when power on, press the button to adjust the speed as you need.

Speed switches every time when press button, shows by 2 LED indicator lights, respectively indicate the Slow and Fast speed.

3. Heating:

It is defaulted to open heating function in startup state.


Press button to close Led light and heating function.

Followed by recycling.

4. Air pressure:

It is defaulted to be in Weak Intensity on power-on state, Press button to switch intensity level to be in Medium, corresponding indicator light on.

1. Power button ON/OFF
2. Heating button: ON/OFF
3. Massage speed (Low / High)
4. Air Compression
(Low / Medium / High)

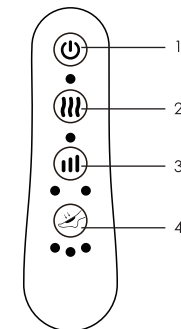	<p>Correct Disposal of this product</p> <p>This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.</p>
--	--

Operations for Use

Problem	Cause	Remedy
Unable to start the device.	The cord is not connected properly.	Please push it harder until you heard the sound of "click".
The appliance stops during use.	The appliance automatically stops after approximately 15 minutes.	This does not indicate a problem.
The main unit does not heat up.	The appliance may not work beyond the recommended temperature.	Use under the recommended temperature of 15 °C - 35 °C

If you still have any problem after using the methods above, please contact your local retailer.

How to clean

1. Cleaning

Before you clean, please make sure the power is switched off, and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.

2. Maintenance

Keep the product clean and far from the wet environment.

Clean when you feel cleaning is necessary.

Cleaning the main unit

- Wipe away dirt with a cloth soaked in soapy water, wring out the cloth before wiping.
- Do not use alcohol, paint thinner, benzine, or nail polish remover. Using these may damage the appliance or cause cracks or discoloration in the components.

Technology parameter

Name: Foot Massager

Model: MaxiGo 2 AE (Series JA)

Rated Voltage: DC12V \approx 2.5A

Rated Power: 30W

DC 24V is not allowed

Manufacturer: XIAMEN EMOKA HEALTHSCIENCE & TECHNOLOGY CO. LTD.

Date of manufacture: look at massager

Made in China

Address of manufacture: D Building No. 100, JINFU ROAD, Chengnan

Industrial Zone, TONG'AN District, Xiamen, China

Phone number: 0086-755-25122982

This product is intended for house-hold use only.

Short time operation 15 minutes, resting time 30 minutes

Note: Certified manufacturer's adapter must be equipped in the use of this products.

General importer in Russia:

OOO "Max Concept", bin 1147453008532

Moscow, Lyubertsy district, Kotelniki, 1st Pokrovsky PR-d, d. 5

Phone number: 8-800-700-14-29

www.ergonova.ru

info@ergonova.ru

Full list of representative offices: <http://ergonova.ru/about/stores/>. All copyrights to texts and images belong to Ergonova GmbH and can not be used otherwise only with the consent of right holder.

Copyright 2012.
Ergonova GmbH
All rights reserved